


Menisoft

Material : Mipafilcon A II 2

DMAA / N-VP / MMA copolymer
(DMAA : Diméthylacrylamide)
72 % water content, non ionic
Aqua-green handling tint
Dk : 34×10^{-11} (cm²/sec) (mlO₂/(ml.mmHg)) ISO 9913-1

Bi-aspheric design

Dk/tc : 42.5×10^{-9} (e_c : 0.08 mm pour -3.00 D) (ISO)

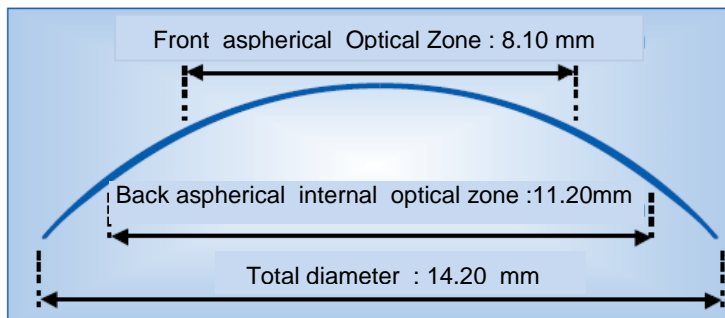

Range

Parameters	r ₀ (mm)	Spherical power (D)	Diameter (mm)
	8.60	-0.25 à -6.00 par 0.25 -6.50 à -10.00 par 0.50	14.20

Fitting rules

One base curve and one diameter, just determine the spherical power.
Wait for 15 minutes after inserting to check the good mobility of the lens.

Indications

- High water content contact lens for myopia correction.
- **Daily wear.**
- **Monthly replacement.**

Lens care

- **MeniCare™Soft**, multipurpose solution for cleaning, disinfection, storage and rinsing.

Packaging

3 lens containers are packaged in a multipack box.
Each Menisoft contact lens is packaged with sterile saline solution in an individual foil sealed tamper-evident plastic lens container, designed to maintain sterility.

